

Contents


STRP-TeraSys - AiO	2
STRP-TeraIMAGE	4
STRP-TeraTune	5
STRP-TeraSys 4000	7
STRP-TeraKit.....	9
Terahertz Generators and Detectors:	11
STRP-KNbO ₃ Crystals	13
Application Notes	14
Publications.....	24

STRP-TeraSys - AiO

The **STRP-TeraSys- AiO** provides a flexible solution for laboratory THz spectroscopy and imaging. It offers maximum flexibility with measurement capabilities in transmission and reflection without realignment of the optics. It is based on organic crystals, to allow access to terahertz frequencies not available with conventional antennas. The **STRP-TeraSys- AiO** includes all optical, mechanical and electronic components for the generation and detection of THz waves such as delay line, terahertz generator, terahertz detector, pump source optics, electronics, humidity sensor, purge chamber, dedicated software, and laptop.



Specifications

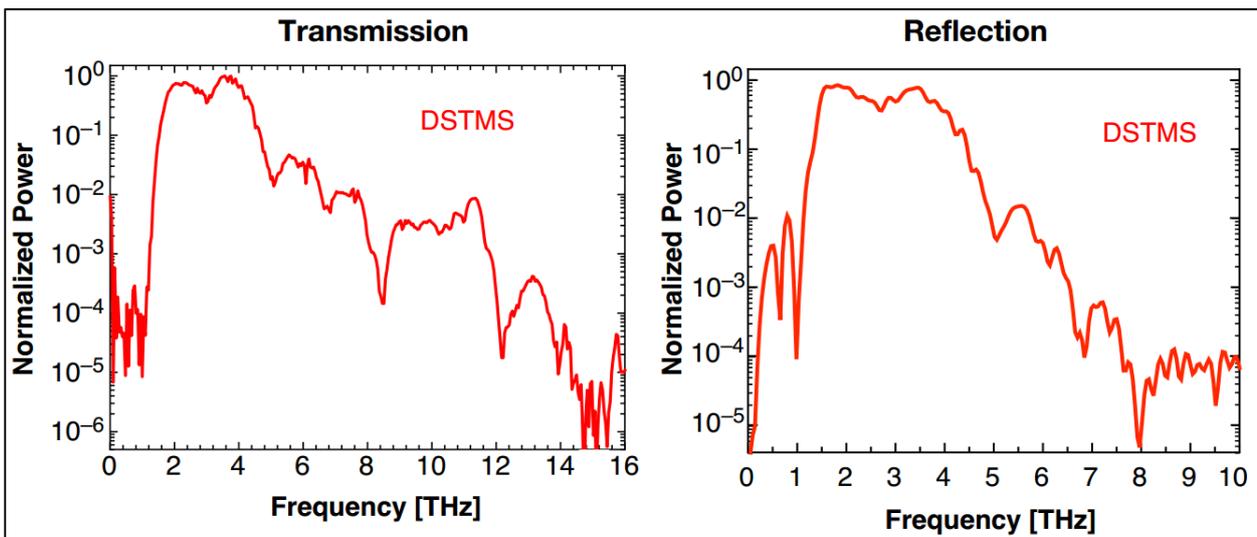

STRP-TeraSys - AiO	Transmission	Reflection
Spectral range	0.3 - 14 THz	0.3 - 8 THz
Dynamic range	> 70 dB	> 40 dB
Signal to Noise (@4 THz)	> 60 dB	> 35 dB
Scan range	up to 60 ps	up to 60 ps
Frequency resolution	< 100 GHz	< 100 GHz
Dimensions (including the pump source)	55 cm x 45 cm x 28 cm	
Pump Source (high power ultrafast Erbium fiber laser)		
Pulse length	< 20 fs	
Total average power	> 200 mW	
Peak power	>120 kW	
Central wavelength	1565 nm	
Repetition rate	80 MHz	

Options

THz imaging with a scanning range of 50 x 50 mm² for transmission for reflection operation. Frequency domain spectrum measured with the **STRP-TeraSys - AiO** using DSTMS as terahertz generator/detector in transmission and reflection.

Applications

Spectroscopy and analysis of materials in the THz range up to 20 THz.

STRP-TeraIMAGE

The flexible solution for THz spectroscopy and imaging

The **STRP-TeraIMAGE** offers a flexible solution for laboratory THz spectroscopy and imaging. It is based on organic crystals, to allow access to THz frequencies not yet available with conventional antennas. The **STRP-TeraIMAGE** includes all optical, mechanical and electronic components for the generation and detection of THz waves such as delay line, THz generator, detector, optics, electronics, dedicated software and laptop. It also features a scanning mechanism for the measurement of phase and full-spectrum images. As the **STRP-TeraKit**, it can also be used with various telecom wavelength lasers.



STRP-TeraIMAGE optical board (Scanning range: 50x50 mm²)

Specifications

THz generator / detector	Organic crystal
Spectral range	1-14 THz
Best phase matchable wavelength	1300-1600 nm

Options

Scanning range of 100x100 mm²


Applications

Spectroscopy and imaging of various materials. Ideal for identification and detection of dangerous substances and biomedical materials, as well as for materials testing of plastic, special polymers and semiconductors. The scan area is 5 x 5 cm². We also offer systems that allow imaging of larger areas.

STRP-TeraTune

The **STRP-TeraTune** provides a tunable narrowband terahertz source with a tuning range from 1.5 to 20 THz and a bandwidth of less than 100 GHz. The basic **STRP-TeraTune** is based on a flash lamp pumped laser (repetition rate 50 Hz to 200 Hz) with a special optical parametric oscillator (OPO) where one wavelength is tuned electronically with a special design to keep the bandwidth below 100 GHz. Other diode pumped laser systems can be offered. The variable wavelength is calibrated and the THz difference frequency is computer controlled (tuning range 1 to 20 THz). Terahertz radiation is generated via frequency mixing in the novel organic nonlinear optical crystal **STRP-DSTMS**, ideally suited to achieve highest terahertz conversion efficiency.

The **STRP-TeraTune** includes all optical, mechanical, electronic components for the generation of THz waves and software control.

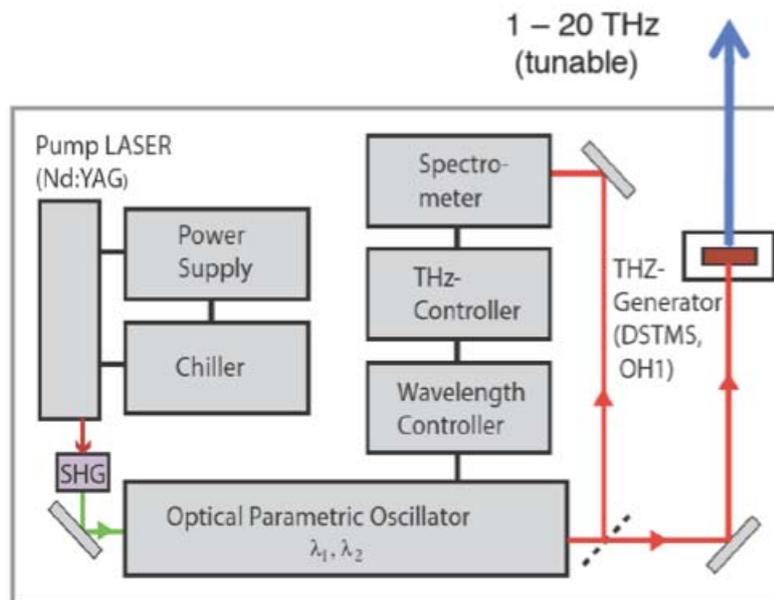


STRP-TeraTune (pump laser, OPO: optical parametric oscillator and THz difference frequency generator, dimensions: 900 x 500 x 180 mm), laser power supply (300 x 450 x 450 mm)

Specifications

- THz frequency tunable 1–20 THz
- Output THz pulse bandwidth < 100 GHz
- Average THz power 1–10 μ W

- THz pulse energy 10–100 nJ
- Linear polarization

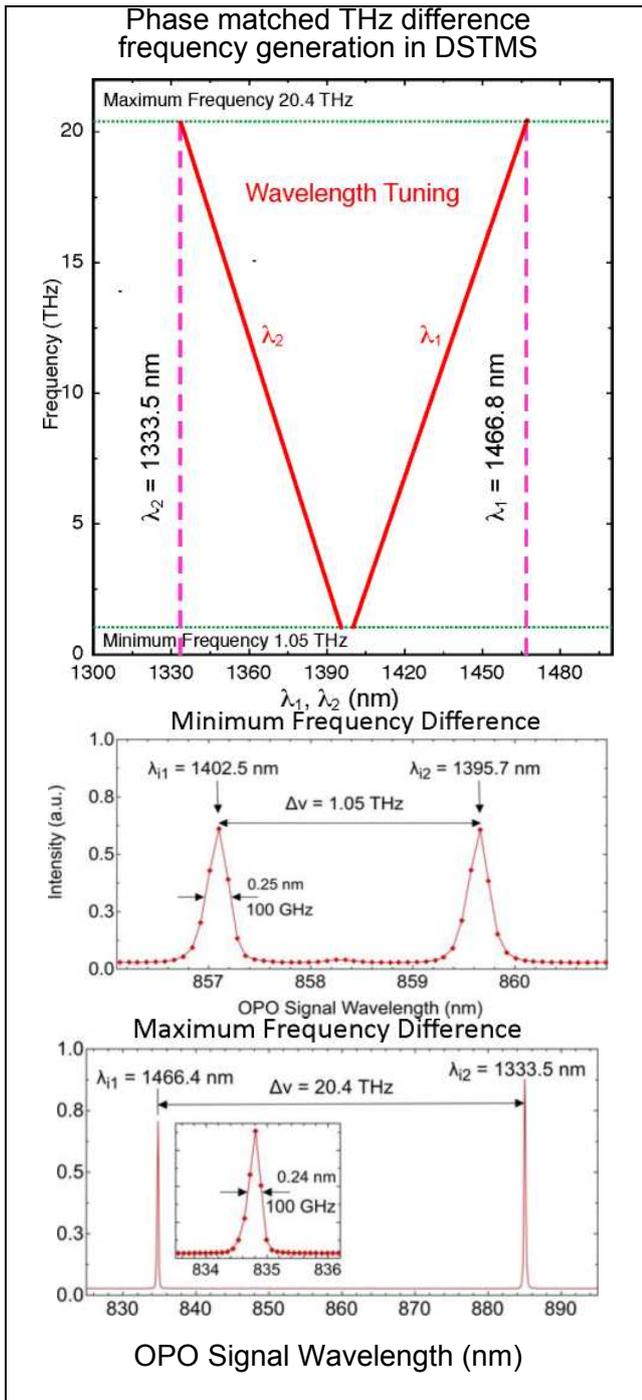


Tunable frequency	
Output frequency	1 - 20 THz
Pulse bandwidth	< 100 GHz
Repetition rate	100 Hz
Average power	1-10 pW
Pulse energy	10-100 nJ
Polarization	linear
OPO: Optical Parametric Oscillator	
Tuning range	1330-1480 nm
Linewidth	< 100 GHz
Repetition rate	100 Hz
Pulse duration	7- 10 ns
Energy per pulse	> 10 mJ
Peak Power	> 1.4 MW
Nd:YAG Pump laser 100-Hz	
Pulse width	10 ns
Pulse energy	> 125 mJ @ 532nm
Divergence	0.5 mrad
Beam quality	$M^2 < 2$
Energy stability	< 1%
Beam diameter	5 mm
Linewidth	90 MHz
Power consumption	< 2.5 kW

STRP-TeraSys 4000

THz System for THz Spectroscopy, Detection and Inspection of Materials

STRP-TeraSys 4000 is a Terahertz spectrometer that operates in the frequency range of 0.3 - 4 THz and has a spectral resolution < 0.01 THz. This turn-key operation system is all solid state, compact and maintenance free. It includes a software package and database for the detection and identification of dangerous materials. It is ideal for applications in spectroscopy, production technology and security.

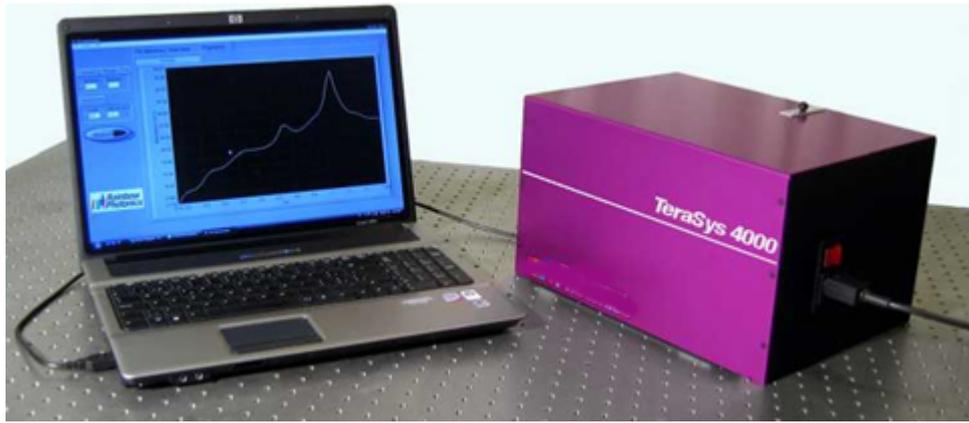

Specification ⁽¹⁾

Frequency range	0.3 - 4 THz
Output power	> 50 nW
Spectral resolution	< 0.01 THz
Polarisation, linear	$> 100: 1$, vertical
Input voltage	110V/ 240V, 50 or 60 Hz
Power consumption	< 60 W
Warm-up time	15 min.
Operating ambient temperature	18°C - 30°C
Dimensions	40x25x18 cm ³


Applications

- Security
 - Explosives detection
 - Bio agents
 - Mail inspection
- Materials inspection
 - Defects in plastics, artificial joints
 - Organic materials
- Spectrography
- Identification of pharmaceuticals, drugs etc.


Features

- Turn key operation

- All solid state
 - Maintenance free
 - Software package and database for the identification of dangerous materials
- ¹ Specifications are subject to change without notice

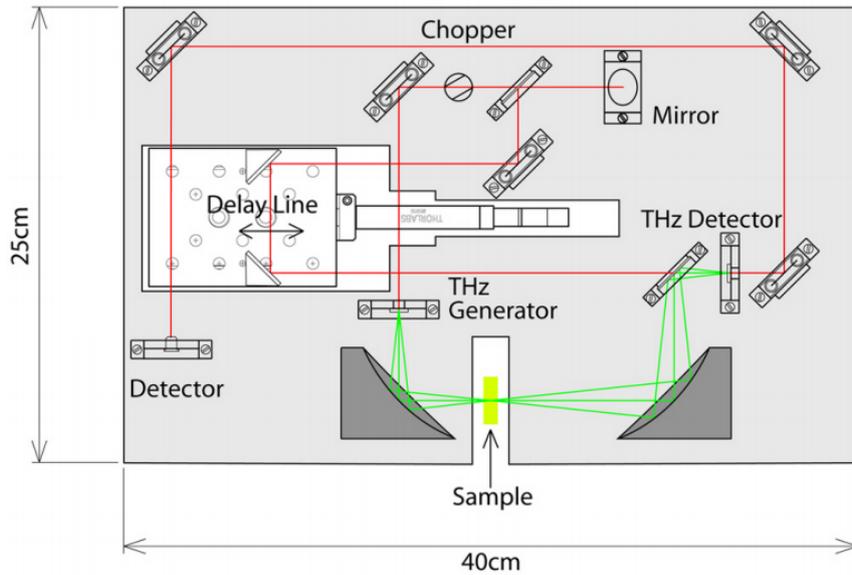

STRP-TeraSys 4000 optical set-up

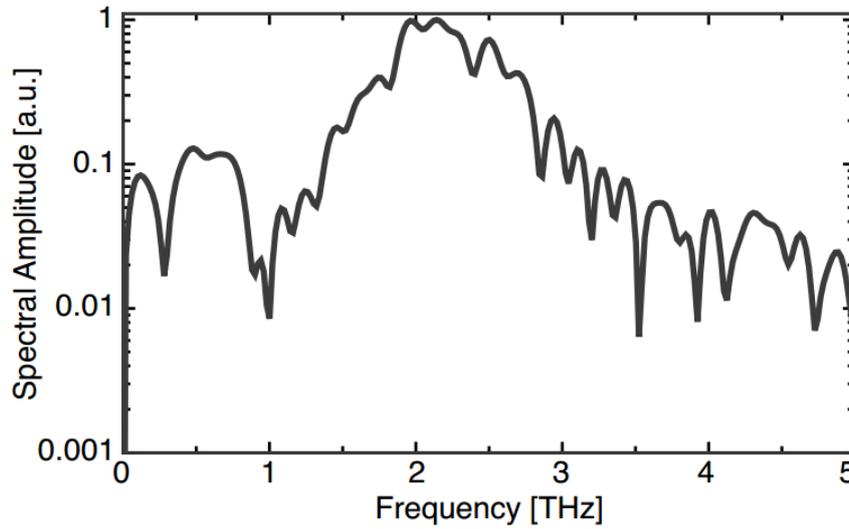

STRP-TeraSys 4000 frequency response.
Other frequency responses are available upon request


STRP-TeraKit

The **STRP-TeraKit** provides a flexible solution for laboratory terahertz spectroscopy. It is based on organic crystals, to allow access to terahertz frequencies by optical down conversion yielding THz frequencies and efficiencies not available with conventional antennas. The **STRP-TeraKit** includes all optical, mechanical, and electronic components for the generation and detection of THz waves such as delay line, terahertz generator, terahertz detector, optics, electronics, lock-in, custom made Er based femtosecond fiber laser, dedicated software and laptop.



STRP-TeraKit optical board (38 cm x 38 cm)

Terahertz spectral bandwidth as a function of the pump pulse length and frequency domain spectrum measured in dry air with the **STRP-TeraKit** using DSTMS as terahertz generator and detector.

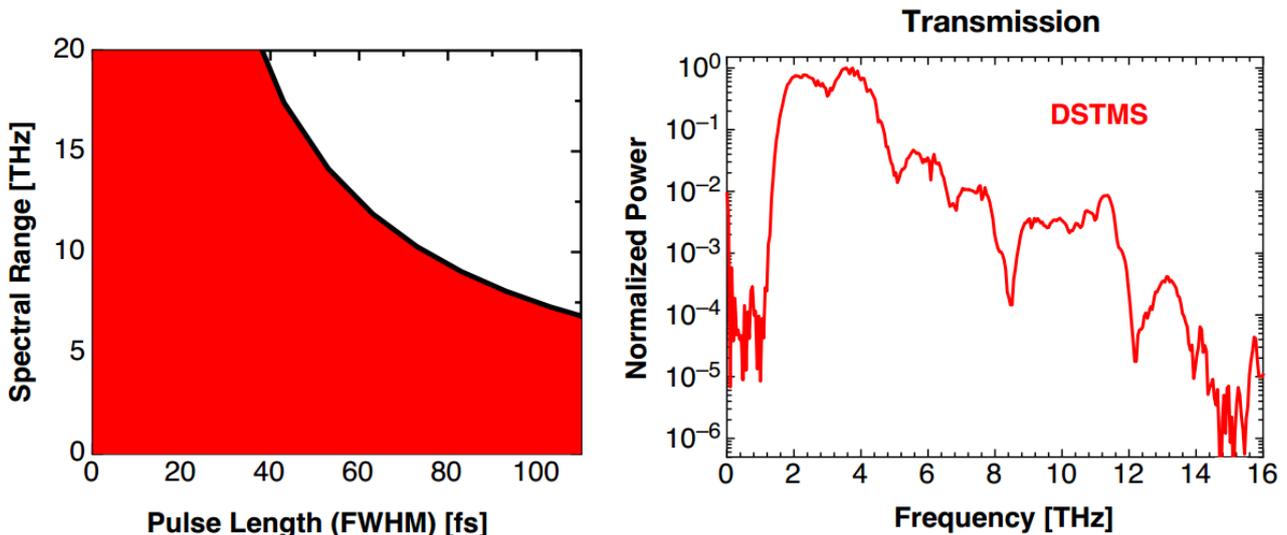

Specifications

THz generator / detector	Organic crystal
Spectral range	1-14 THz
Best phase matchable wavelength	1300-1600 nm
Dynamic range	> 60 dB, (@4THz)
Scan range	up to 60 ps

Frequency resolution	< 100 GHz
Dimensions	30 x 38 x 17 cm
<u>Pump source (Er based femtosecond laser)</u>	
Pulse length	< 20 fs
Total average power	> 200 mW
Peak power	> 120 kW
Central wavelength	1565 nm
Repetition rate	> 80 MHz

Options

THz imaging with a scanning range of 50x50 mm² or 100x100mm²

Applications

Spectroscopy and analysis of materials in the THz range up to 20 THz (depending on the pump laser).

Other spectral ranges are available upon request.

Terahertz Generators and Detectors:

STRP-OH1 Crystals: (2-(3-(4-Hydroxystyryl)-5, 5-Dimethylcyclohex-2-Enylidene)Malononitrile)

THz generators/detectors are based on various organic crystals: DAST, DSTMS and OH1. They are optimized for operation at Terahertz frequencies from 0.3 - 20 THz. The organic electro-optical crystals (DAST, DSTMS, OH1) have high electro-optic coefficients and show ultra-fast electro-optical response for modulation frequencies up to 200 GHz. These crystals can also be used for electro-optic detection.



Example of an as-grown **STRP-OH1** crystal and organic crystals polished, coated and mounted for optical applications. **STRP-OH1** and other high-quality organic crystals are produced and optically prepared in the facilities of Rainbow Photonics in Switzerland.

Properties

- High quality crystals
- Cut and polished for various applications
- Large nonlinear optical susceptibilities
- Large electro-optic coefficients
- Phase matching for THz-wave generation between 1200 nm and 1460 nm

Features

- High quality crystals
- Cut and polished for various applications
- Large nonlinear optical susceptibilities
- Phase matching for THz generation between 1200 nm and 1460 nm


Applications

- Efficient THz generation and detection from 0.1 to >10 THz
- Fast electro-optic modulation
- Optical parametric generation
- Efficient frequency doubling of 1.55 μm radiation


Physical Properties

Melting point	212 °C
Point group symmetry	mm2
Refractive indices	$n_2 = 1.58, n_3 = 2.15$
Nonlinear coefficients ($A=1.9 \text{ } \mu\text{m}$)	$d_{333} = 120 \pm 10 \text{ pm/V}$ $d = 13 \pm 2 \text{ pm/V}$ $d_{322} = 8.5 \pm 2 \text{ pm/V}$
Electro optic coefficients	$r_{333} (633\text{nm}) = 109 \pm 4\text{pm/V}$ $r_{333} (785 \text{ nm}) = 75 \pm 7\text{pm/V}$ $r_{333} (1064\text{nm}) = 56 \pm 2\text{pm/V}$ $r_{333} (1319\text{nm}) = 52 \pm 7\text{pm/V}$

Absorption Spectrum

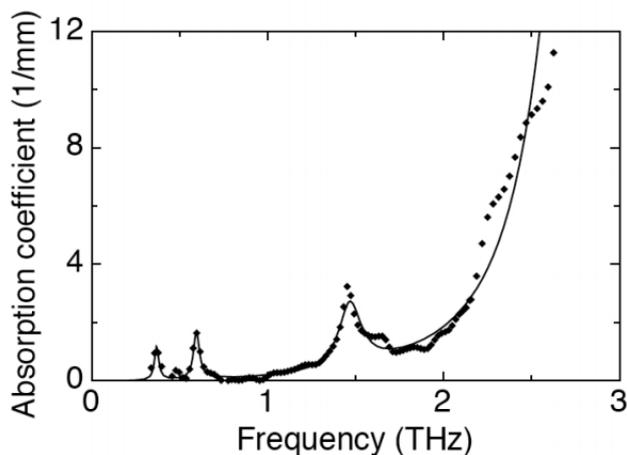

THz Conversion Efficiency

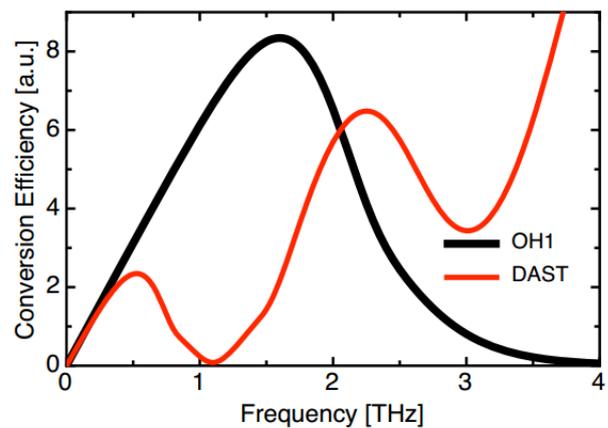


- 1) "Configurationally locked, phenolic polyene organic crystal OH1: linear and nonlinear optical properties"; C. Hunziker, S. Kwon, H. Figi, F. Juvalta, O. Kwon, M. Jazbinsek, P. Gunter, J. Opt. Soc. Am. B 5, 1678 (2008).
- 2) "A hydrogen-bonded organic nonlinear optical crystal for high-efficiency terahertz generation and detection"; F. Brunner, O. Kwon, S. Kwon, M. Jazbinsek, A. Schneider, P. Gunter, Opt. Express 16, 16496 (2008).
- 3) "Organic phenolic configurationally locked polyene single crystals for electro-optic and terahertz wave applications"; O. Kwon, S. Kwon, M. Jazbinsek, F. Brunner, J. Seo, C. Hunziker, A. Schneider, H. Yun, Y. Lee, P. Gunter, Adv. Funct. Mater. 18, 3242 (2008).
- 4) "Scaling submillimeter single-cycle transients toward megavolts per centimeter field strength via optical rectification in the organic crystal OH1"; C. Ruchert, C. Vicario, C.P. Hauri, Opt. Lett. 37, 899 (2012).

STRP-KNbO₃ Crystals

Standard Features

- High quality nominally undoped crystals for electro-optics and nonlinear optics
- High quality crystals doped with Rh, Fe, Mn, and Ni for photorefractive applications in the visible and infrared
- Very low scattering losses

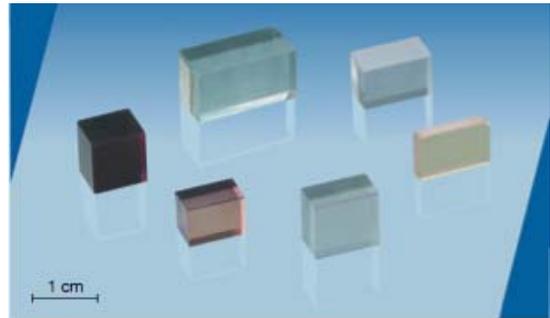

Optional Features

- High photosensitivity up to 1000 nm wavelength
- Milliseconds response time

Applications

- Electro-optics and nonlinear optics
- Photorefractive applications with laser diodes


Dynamic holography and optical phase conjugation in the visible and near infrared

Photorefractive Grating Recording Times

(Selected STRP- KNbO₃ crystals at different wavelengths for $I = 1 \text{ W/cm}^2$)

	Wavelength [nm]	Recording time (typical) [s]
STRP-KNbO ₃ :Fe	488	1
STRP-KNbO ₃ :Mn	515	1
	860	3
STRP-KNbO ₃ :Fe reduced	488	0.01
	515	0.01
STRP-KNbO ₃ :Rh reduced	860	0.5
	1064	50

Absorption Spectra

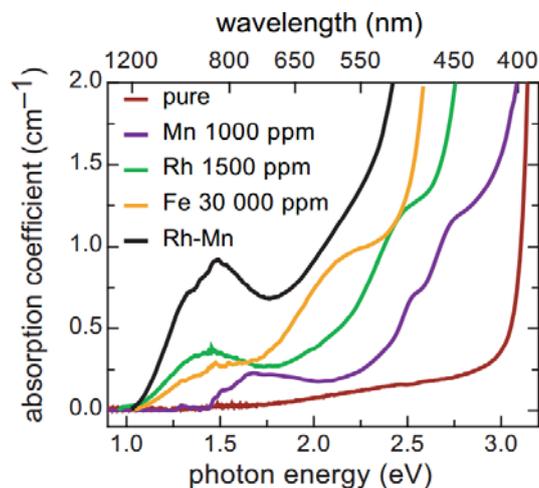

Application Notes

Broadband terahertz generation and detection Widely tunable and narrowband THz sources Applications of THz waves

The interest in terahertz electromagnetic radiation stems from the unique interactions of these rays with matter, which can be exploited in various applications. Terahertz waves excite molecular vibrations and lattice vibrations in materials, which make THz radiation very interesting for spectroscopy and material identification.

Terahertz radiation is non-ionizing, very sensitive to water and hydration state, and transparent to non-polar substances such as nonconductive polymers, paper, packing material, etc. Therefore material irregularities, defects, enclosures, not seen with other techniques, can be visualized with terahertz radiation. THz radiation has therefore a high potential for the non-destructive materials testing.

There are several techniques to produce terahertz radiation. We have developed **novel organic THz generators**, with optimized properties that make them ideal materials for the efficient generation and detection of terahertz radiation. Based on these materials we developed compact systems for THz time-domain spectroscopy, **STRP-TeraSys** and **STRP-TeraKit**, as well as for THz imaging, **TeraIMAGE**. We have also developed a unique tunable single-frequency THz source, **TeraTune**, with a very wide tunability range 120 THz and a narrow linewidth < 100 GHz.

Terahertz radiation

The terahertz range of the electromagnetic spectrum is situated between high frequency electronics (microwaves) and long wavelength photonics (infrared light). Terahertz radiation is easily accessible by means of blackbody radiation, but it is a challenge to separate the signals from the natural background. Several techniques have been used to generate electromagnetic waves below and above this so-called "THz frequency gap". Figure 1 shows the rough power levels of some of them. Electronic techniques can be used to generate waves with frequencies (mainly by electronic frequency multiplication of lower frequency sources) up to about 0.5 THz. From 0.3 to 3 THz Auston switches are very popular sources. Nonlinear optical techniques (optical rectification and difference-frequency generation) can be used to cover the frequency range between 0.3 to 50 THz and quantum-cascade lasers between roughly 20 THz and 100 THz.

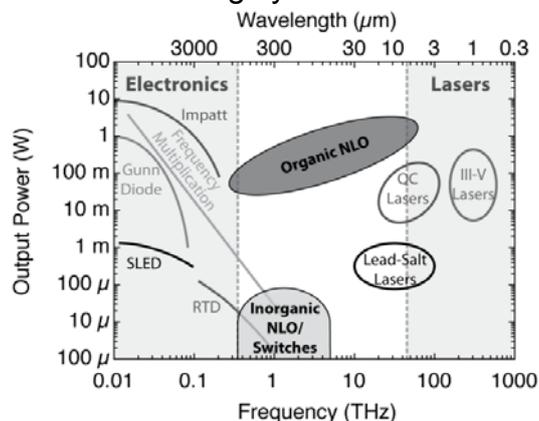


Figure 1: Approximate output power (either CW or peak power) as a function of frequency (wavelength) of electromagnetic radiation in the THz range and around when using different generation techniques. Adapted and updated from M. Tonouchi, Nat. Photonics 1, 97 (2007).

Broadband THz sources

Most broadband THz sources are based on the excitation of different materials with ultrashort laser pulses in the femtosecond range. Photoconduction and optical rectification are two of the most common approaches for generating broadband THz pulses. Optical methods are used for the generation of broadband terahertz sources, and due to the increasing progress in laser technology these methods have been the most developed in the last 20 years.

In the **photoconductive** approach a femtosecond laser generates an ultrafast photocurrent in a photoconductive switch or semiconductor using electric-field carrier acceleration. The achievable bandwidth is limited to a few THz due intrinsic limits of the carrier lifetime in semiconductors.

Optical Rectification is an alternative mechanism for pulsed THz generation. A femtosecond laser is used as well, and the energy of the terahertz radiation comes directly from the exciting laser pulse. In this case the conversion efficiency depends, beside the parameters of the pump laser, mainly on the electro-optic coefficient and the proximity of the phase-matching conditions of the material.

In optical rectification, a high-intensity ultrashort laser pulse passes through a transparent crystal material that emits a terahertz pulse without any applied voltages. Figure 2 shows a schematic of optical rectification with a pulsed femtosecond laser with organic crystal generators DAST or DSTMS. In this nonlinear-optical process, a nonlinear material is quickly electrically polarized at high optical intensities. This changing electrical polarization emits terahertz radiation. It is called rectification because the rapid oscillations of the electric field of the laser pulse are "rectified" and only the envelope of the oscillations remains. Since the medium absorption is low, the polarization instantaneously follows the pulse envelope implying that there is practically no limit on the speed at which the polarization can be switched on and off, i.e. there is no intrinsic limit on the bandwidth as for the photoconductive antennas.

The bandwidth of the THz pulse generated by the optical rectification process is limited by the length of the pump laser and the velocity matching conditions in the generator and detector material. THz bandwidth achievable for pump laser sources with various pulse lengths is illustrated in Figure 3. The materials for THz generation by optical rectification are discussed in the following section.

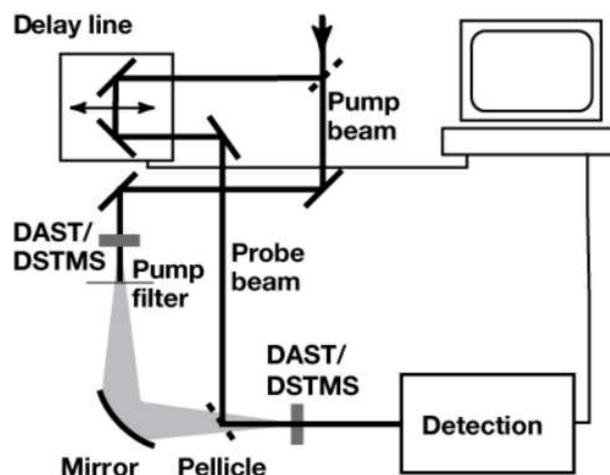


Figure 2: Scheme of the setup for generation and detection of THz pulses in organic THz generators DAST or DSTMS. The generation is based on optical rectification of femtosecond pulses, while the detection is based on electro-optically induced

changes of the refractive index in a THz detector, induced by the THz electric field and probed by an optical probe beam. For inorganic semiconductors most often the electro-optic sampling is used,¹ while we use THz induced lensing in case of organic crystals.²

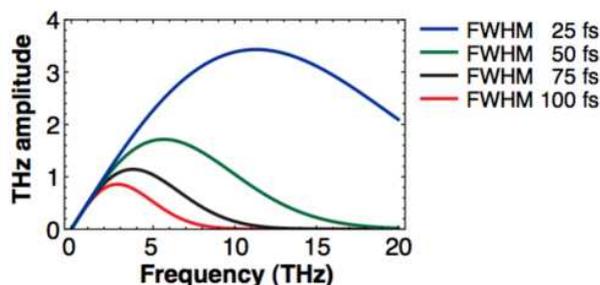


Figure 3: Terahertz frequency range for femtosecond laser sources with various pulse lengths (FWHM from 25 fs to 100 fs) in case of perfectly phase-matched THz generation materials.

Materials for THz generation

Due to the larger nonlinear optical susceptibilities and velocity matching between THz and optical pump waves of organic materials compared to the inorganic ones, much larger power levels, limited by the damage thresholds of the materials, can be obtained by using organic materials as THz generators.

Table 1 shows most relevant material parameters for THz wave generation for best organic crystals compared to inorganic crystals as well as for an electro-optic polymer LAPC. As it can be seen in this table, the organic crystals OH1, DSTMS and OH1 show the largest figure of merit FM_{THz} . Phase matching is possible when the THz refractive index n_{THz} is close to the group index n_g at the mean optical pump wavelength.³ This kind of phase matching is usually referred to as group velocity matching and is commonly used for broadband THz-wave generation by optical rectification of fs pump lasers.³

	n_o	n_g	n_{THz}	r (pm/V)	FM_{thz} (pm/V) ²	Table 1: Relevant parameters of organic and inorganic nonlinear optical materials that have been investigated for optical-to-THz frequency conversion. Refractive index n_o and the group index n_g at the pump optical wavelength; refractive index n_{THz} in the THz range; electro-optic coefficient r ; figure of merit FM_{THz} for THz generation by optical rectification. Organic crystals: DAST, DSTMS, OH1 LAPC: guest-host polymer [X.M. Zheng et al, J. Nanoelectron. 2, 58 (2007)] Inorganic crystals: GaAs, ZnTe, InP, GaP, ZnS, CdTe, LiNbO3
DAST	2.13	2.26	2.29	47	5600	
DSTMS	2.08	2.19	2.20	49	6100	
OH1	2.16	2.33	2.28	52	7400	
LAPC	1.6	1.8	1.7	52	1700	
GaAs	3.37	3.61	3.63	1.6	66	
ZnTe	2.83	2.18	3.16	4	160	
InP	3.2	3.16	3.54	1.45	40	
GaP	3.12		3.34	1	17	
ZnS	2.3		2.88	1.5	7	
CdTe	2.82		3.24	6.8	470	
LiNbO ₃	2.2	2.18	4.96	28	1100	

$$FM_{THz} = \frac{n_o^6 r^2}{16n_{THz}}$$

¹ G. Gallot and D. Grischkowsky, "Electro-optic detection of terahertz radiation," J. Opt. Soc. Am. B 16, 1204 (1999)

² A. Schneider, I. Biaggio, P. Gunter, "Terahertz-induced lensing and its use for the detection of terahertz pulses in a birefringent crystal," Appl. Phys. Lett. 84, 2229 (2004)

³ Schneider, A.; Neis, M.; Stillhart, M.; Ruiz, B.; Khan, R. U. A. & Gunter, P., "Generation of terahertz pulses through optical rectification in organic DAST crystals: theory and experiment," J. Opt. Soc. Am. B, 2006, 23, 18221835

Best inorganic electro-optic materials, such as LiNbO_3 are far from optimal phasematching conditions and can only be used in special configurations. Therefore, although having lower electro-optic coefficients and figures of merit, most commonly semiconducting electro-optic materials, such as ZnTe are used because they can operate close to phase matching. Organic materials combine both high figures of merit and possibility for velocity matching, therefore we chose these materials for our instruments. Organic crystals DAST, DSTMS and OH1 show best phase matching for fs laser sources in the telecom wavelength range 1200-1600 nm, which makes them very attractive for compact table-top THz instruments. They can also operate close to velocity matching in a broad range of THz frequencies, making possible broadband THz generation with low-power fs laser sources. Figure 4 shows an example of a broadband THz field generated in organic crystal DSTMS using **STRP-TeraKit** in comparison with the field generated by a semiconducting antenna using a compact fs fiber laser source.

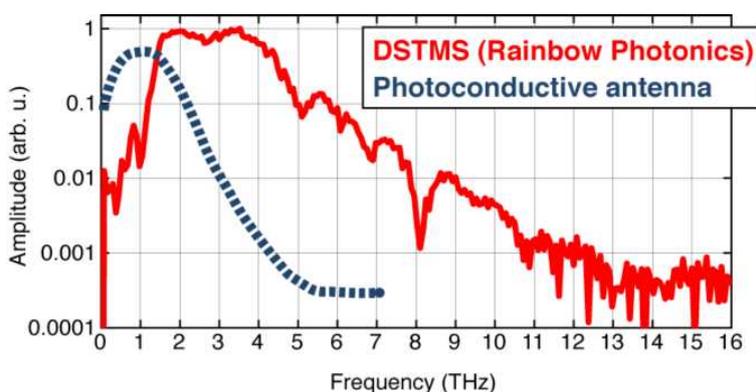


Figure 4: THz field amplitude as a function of frequency for a THz pulse generated in an organic crystal DSTMS (in STRP-TeraKit) using femtosecond pump lasers and THz time-domain spectroscopy, and typical range reached with PC antennas.

Organic THz generators and detectors



Figure 5: THz generators/detectors optically prepared and mounted for applications. Standard aperture sizes range from 2 to 5 mm; larger sizes are available upon request.

The generators exhibit a high damage threshold of over 150 GW/cm^2 at 150 fs pulse length and 1500 nm.

High THz electric fields of over MV/cm have been generated by our generators.^{4,5}

Terahertz generators and detectors are produced and optically prepared at the facilities of our Photonics. We are the worldwide only commercial producer of organic single crystalline THz generators.⁴

⁴ Hauri, C. P.; Ruchert, C.; Vicario, C. & Ardana, F., "Strong-field single-cycle THz pulses generated in an organic crystal," Appl. Phys. Lett. 99, 161116 (2011)

⁵ Ruchert, C.; Vicario, C.; Hauri, C. P., "Scaling submillimeter single-cycle transients toward megavolts per centimeter field strength via optical rectification in the organic crystal OH1," Opt. Letters 37, 899 (2012)

Figure 6 shows best application ranges of our THz generators: pump lasers in the wavelength range 1200-1600 nm are most suitable. Shorter the pulse, larger will be the generated THz bandwidth when using optical rectification, as shown in Figure 3.

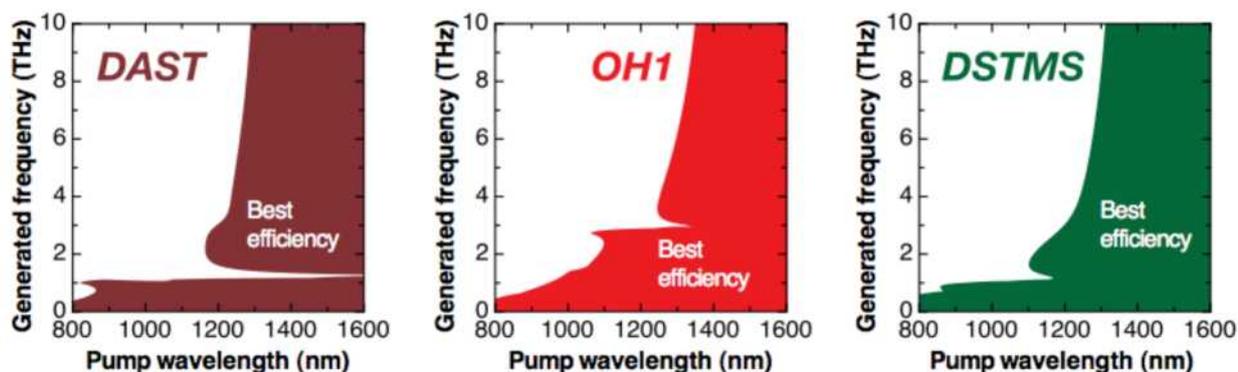


Figure 6: Best application ranges of organic single-crystalline THz generators DAST, OH1 and DSTMS. The shaded ranges present areas of coherence length exceeding 0.5 mm. Pulsed lasers with wavelengths between 1200 and 1600 nm are best suited as pump lasers for our THz generators.

Detection: In THz time domain spectroscopy, one can use common techniques for the detection of THz signals generated in organic crystals. To achieve best detection efficiency in a broad THz range, same organic materials can be also used for electrooptic detection using a different principle as for standard electro-optic sampling,¹ which is limited to optically isotropic materials. For organic crystals (or other birefringent materials) THz-induced lensing can be used with a similar sensitivity as the electro-optic sampling.²

Terahertz time-domain spectroscopy and imaging with *STRP-TeraKit* and *STRP-TeraIMAGE*

The **STRP-TeraKit** provides a flexible solution for laboratory terahertz spectroscopy. It is based on organic crystals, to allow access to terahertz frequencies not yet available with conventional antennas. The **STRP-TeraKit** includes all optical, mechanical and electronic components for the generation and detection of THz waves such as delay line, terahertz generator, terahertz detector, optics, electronics, dedicated software and laptop. It can be used with any femtosecond laser source at telecom wavelengths.

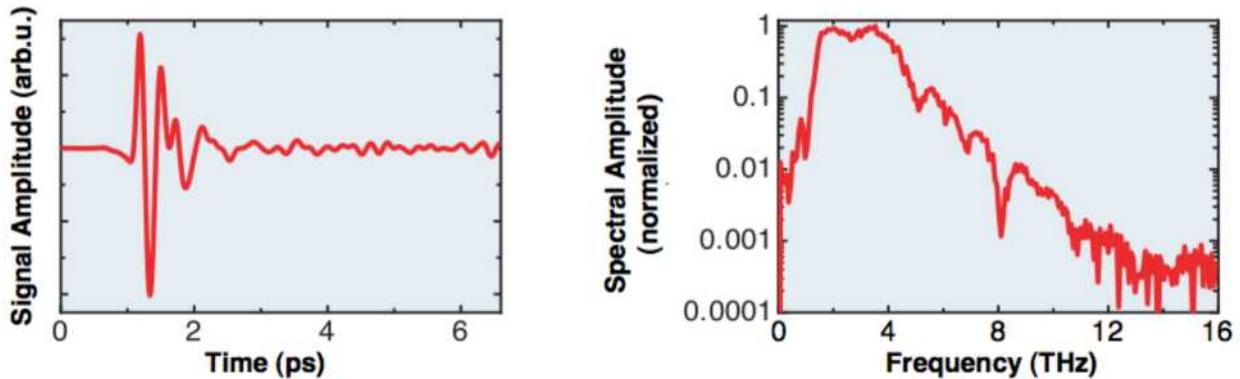


Figure 7: *STRP-TeraKit* optical board (38cm x 38cm).




Figure 8: *TeraIMAGE* optical board (scanning range 50mm x 50mm).

The THz time-domain spectrometer with imaging option, *STRP-TeraIMAGE*, includes, Figure 9: THz time-domain signal amplitude and amplitude spectrum in *STRP-TeraKit*, using 0.45 mm thick The THz detection in *STRP-TeraKit* is optimized using special optical and electronic components (the details are confidential) so that a high SNR (signal to noise ratio) can be achieved already using relatively low-power femtosecond fiber lasers. An example of the time-domain THz signal and its spectrum is shown in Figure 9.



DSTMS crystals for THz generation and detection, and a pump laser with 65 fs pulse length, average power of 180 mW and 1.8 nJ energy/pulse.

Optical image:

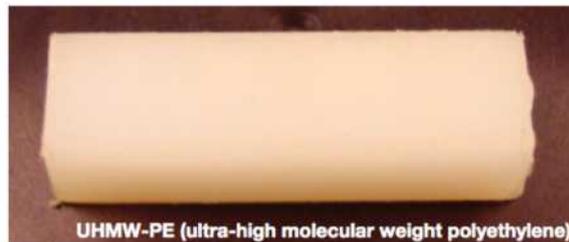

THz image:

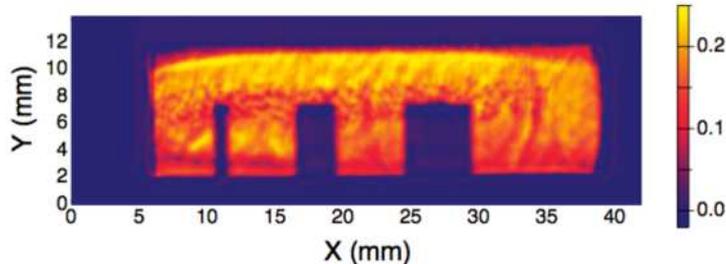


Figure 10: Optical image (made by a usual photocamera) and terahertz image (made by *STRP-TeraIMAGE*) of a piece of plastics with optically invisible defects.

The THz time-domain spectrometer with imaging option, *STRP-TeraIMAGE*, includes, beside the spectroscopic part, which is the same as in *STRP-TeraKit*, the imaging part with all necessary mechanics control and data acquisition software for scanning objects up to 50 mm x 50 mm (larger available upon request). An example of THz imaging with *STRP-TeraIMAGE* is shown in Figure 10.

Narrowband widely tunable THz source: *STRP-TeraTune*

Many materials exhibit specific absorption features (fingerprints) not only in the THz range up to about 3 THz, mostly investigated because this range that can be reached by photoconductive antennas, but also higher above, see Table 2. Also, the attenuation in air due to water-vapor absorption, which limits the range of application possibilities below about 10 THz becomes much smaller (up to four orders of magnitude at 18 THz!) above 10 THz (see Figure 11), which makes it interesting to extent the THz range up to 20 THz.

Table 2: THz fingerprints of explosives and related compounds [from Liu et al., Proc. IEEE].⁶

Explosive & Related Compound	Measured Absorption Peak Position (THz)
TNT	1.66, 2.20, 3.69, 4.71, 5.52, 8.28, 9.12, 9.78, 10.65, 11.01, 13.86, 15.15, 16.95, 17.37, 19.17, 19.89
RDX	0.82, 1.05, 1.50, 1.96, 2.20, 3.08, 6.73, 10.35, 11.34, 12.33, 13.86, 14.52, 17.74, 18.12, 20.13
HMX	1.78, 2.51, 2.82, 5.31, 6.06, 11.28, 12.00, 12.54, 12.96, 13.74, 14.55, 18.15, 18.60, 19.38
PETN	2.0, 2.84
Tetryl	5.97, 10.11, 11.28, 14.67, 16.14, 18.36
2-amino-4, 6-DNT	0.96, 1.43, 1.87, 3.96, 5.07, 6.27, 8.49, 9.87, 10.77, 12.15, 13.44, 16.68
4-amino-2, 6-DNT	0.52, 1.24, 2.64, 3.96, 5.04, 5.82, 7.53, 9.30, 10.20, 11.13, 13.86, 14.97, 17.70
4-Nitrotoluene	1.20, 1.37, 1.86, 6.75, 8.85, 10.83, 14.04, 15.66, 18.51
1,3,5-TNB	4.17, 4.62, 10.05, 11.19, 13.80, 15.75, 19.05
1,3-DNB	0.94, 1.19, 2.37, 10.56, 12.18, 15.33, 17.13
1,4-DNB	3.24, 3.5, 5.55, 10.38, 12.45, 13.29, 15-21, 15.54
2,4-DNT	0.45, 0.66, 1.08, 2.52, 4.98, 8.88, 10.56, 11.58, 12.81, 14.34, 15.69, 19.05, 20.04
2,6-DNT	1.10, 1.35, 1.56, 2.50, 5.61, 6.75, 9.78, 11.43, 13.32, 13.89, 15.39, 17.25
3,5-dinitro aniline	0.96, 1.20, 3.18, 4.62, 5.04, 5.91, 7.44, 10.62, 10.98, 14.46, 16.41, 18.18
2-nitro diphenyl aniline	2.19, 2.58, 2.88, 3.45, 5.13, 6.18, 7.56, 10.08, 12.33, 13.05, 15.00, 15.60, 16.29, 17.34, 18.51, 19.32

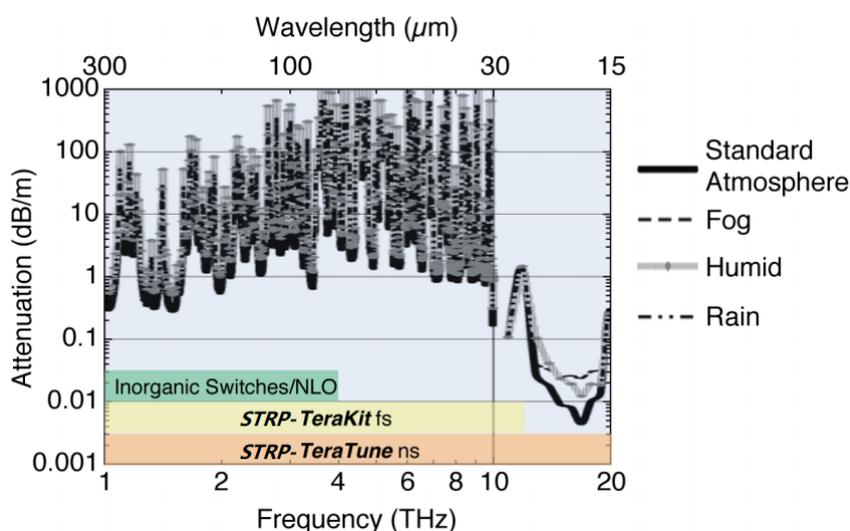


Figure 11: Attenuation of THz radiation in air [from Appleby et al, IEEE 2007]⁷⁷. Above 10 THz, the attenuation in air dramatically decreases compared to the attenuation below 10 THz. This opens up a wide range of new application possibilities, such as remote imaging and sensing.

For some applications a high THz beam power in a narrow band can be preferable than a broadband pulse. The total THz power generated with broadband generation techniques is distributed over the spectral content of the pulses; therefore the power density at any particular frequency is inherently low. To obtain reasonable conversion efficiency for a certain THz frequency, a narrower band pulsed output with high beam peak power is preferred.

⁶ H.B. Liu, H. Zhong, N. Karpowicz, Y. Chen, X.C. Zhang, "Terahertz Spectroscopy and Imaging for Defense and Security Applications," Proc. IEEE 95, 1514 (2007).

⁷ R. Appleby, H. B. Wallace, "Standoff detection of weapons and contraband in the 100 GHz to 1 THz region", IEEE Trans. Antennas Propag. 55, 2944 (2007)

STRP-TeraTune is a tunable narrowband THz source with a narrow linewidth of less than 100 GHz and a large tuning range 1-20 THz. We introduced this unique THz source to the market in 2012. It is based on difference frequency generation of nanosecond pulses in organic THz generators DSTMS and OH1. The suitable infrared pump pulses are generated in a specially designed dual-wavelength OPO (optical parametric oscillator), which is tunable in the range of 1330-1480 nm and produces two narrowband nanosecond pulses with their frequency difference in the THz range. The wavelengths can be tuned by angular tuning of OPO crystals, which is motorized and controlled by the accompanying software.

The generator crystals can be externally switched for achieving the best efficiency at different THz frequencies. High THz peak power of more than 30 W can be reached at 1.25 THz using OH1 generator crystals. The tunability spectrum using 1-mm thick DSTMS crystals is shown in Figure 13.

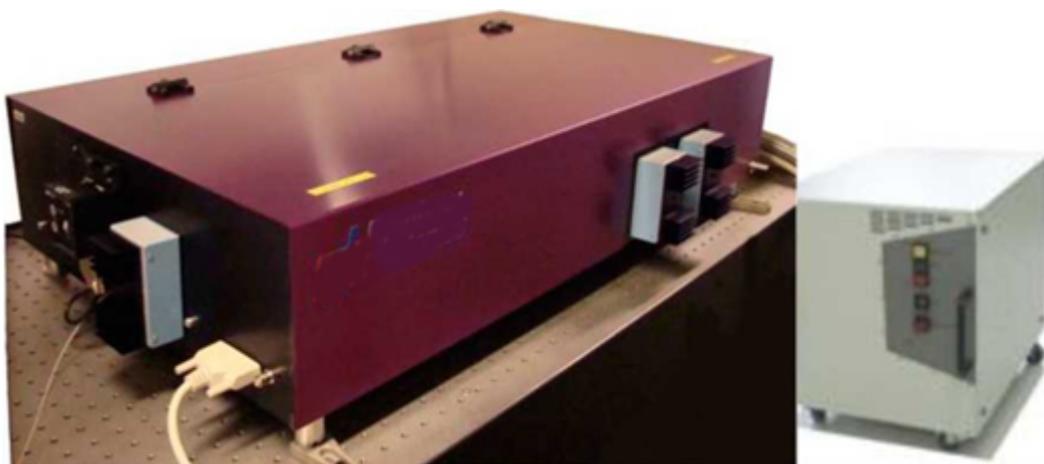


Figure 12: *STRP-TeraTune*: Tunable (1-20 THz) narrow linewidth (<100 GHz) THz source.

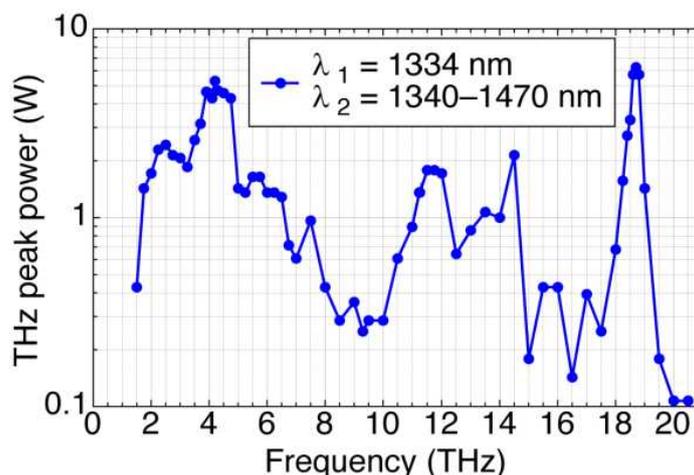


Figure 13: Tuning curve of *STRP-TeraTune* using a 1-mm DSTMS THz generator. In some regions the achievable THz peak power is low, which is either due to the THz absorption of the generation material itself or due to non-perfect phase matching. The generation crystal is relatively thick (1 mm), so that high efficiency is achieved where the generation is phase matched.





Applications of THz waves

Some of the applications of THz waves are related to the unique properties of these waves to excite molecular vibrations and lattice vibrations in the "Reststrahlen" range. In addition the THz waves show low absorption and are transmitted through most non-conductive

homogeneous plastics, paper, carton, most clothes, etc, and can therefore detect hidden hazardous substances. Therefore besides THz spectroscopy of materials these waves are potentially useful for security applications, but also for the identification of defects in non-conductive materials. For conductive and partially conductive materials, THz spectroscopy can give useful insights into the mechanisms of charge transport in these materials. Here we give some examples of THz spectra and materials testing demonstrated using organic nonlinear optical materials for generation and detection of THz waves.

Figure 14(a) shows the THz spectra of several explosives as measured using THz time-domain spectroscopy.

Figure 14(b) shows a Semtex explosive sample hidden behind two teflon plates as seen by optical waves (left) and THz waves (right; Semtex: green; the yellow area on the upper right corresponds to the red paper sticker shown in the left picture).

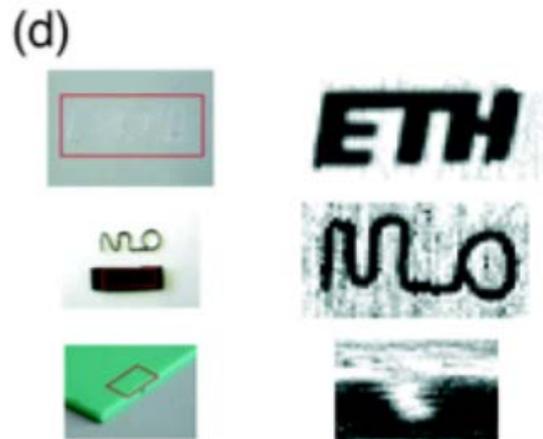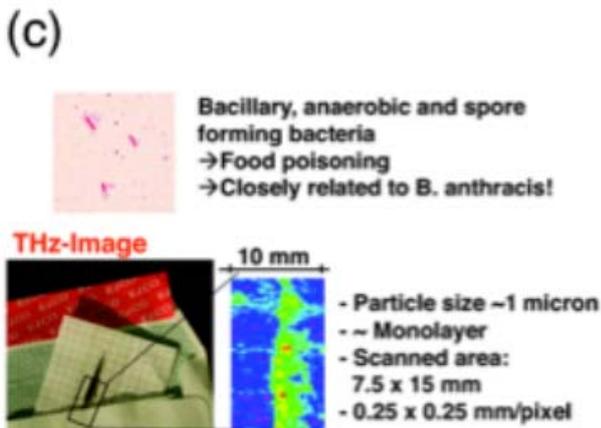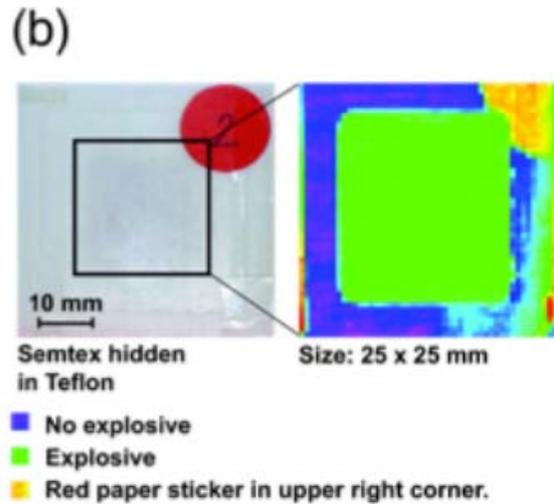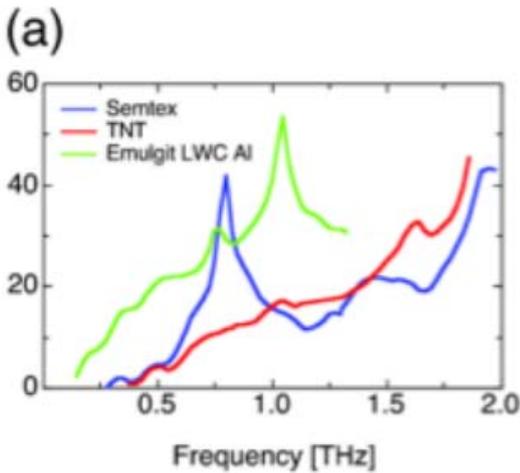

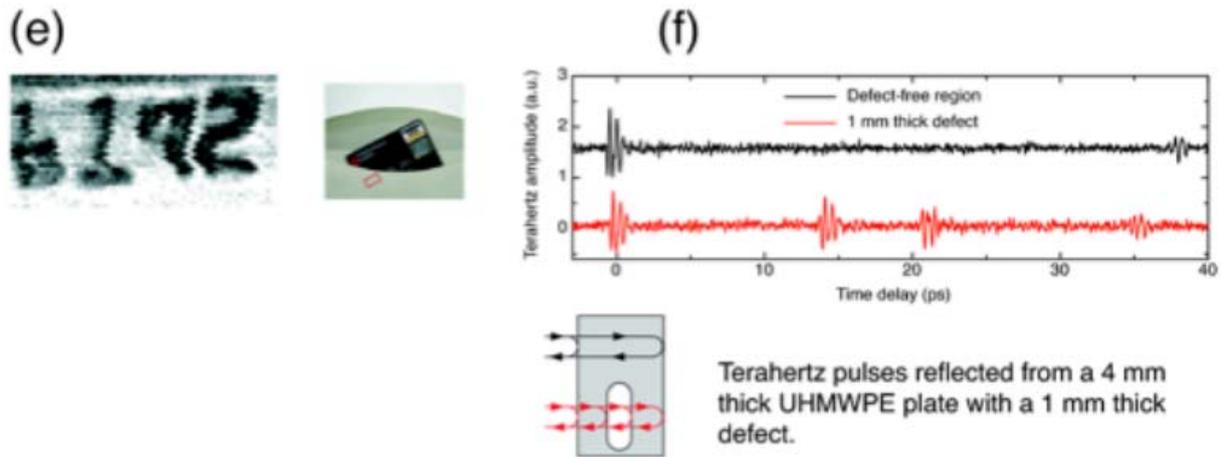


Figure 14: Examples of applications demonstrated by using organic nonlinear optical materials for THz generation and detection.

Figure 14(c) shows a picture and THz image of a bacillus cereus spores (anthrax) sample hidden in an envelope. Examples of materials testing are shown in Figure 14(d) and (e). Figure 14(d) shows the optical and THz pictures of a pile of overhead transparencies, where the label "ETH" has been cut out in one of the transparencies (not seen by visible light), and its THz transmission image giving a full contrast image due to the phase shift of the THz wave in the cut-out area (the "defect"). The second picture shows how metallic defects or inclusions (a metal wire with the symbol "NLO") embedded in a plastic can be made visible by THz waves.

The last picture shows the identification of a void in a piece of plastics and Figure 14(e) the impressed credit card number from a credit card hidden in an envelope.

Figure 14(f) shows examples of THz reflection images of polyethylene samples with and without defects. The voids can be made visible in 3D with a resolution of less than 10 pm. This lower than the wavelength (longitudinal) resolution is due to the fact that the phase shift and temporal resolution of the reflected wave can be determined very precisely.

Publications

THz generation and detection using terahertz instruments

1. Andreja Abina, Uros Puc, Anton Jeglic, Aleksander Zidansek
"Structural characterization of thermal building insulation materials using terahertz spectroscopy and terahertz pulsed imaging" Non-Destructive Testing and Evaluation International 2016, 77, 11-18
2. Majkic et al, "Optical properties of aluminum nitride single crystals in the THz region", Optical Materials Express 5, No. 10 (2106) 2015
3. Uros Puc, Andreja Abina, Melita Rutar, Aleksander Zidansek, Anton Jeglic, Gintaras Valusis "Terahertz spectroscopic identification of explosive and drug simulants concealed by various hiding techniques" Applied Opt. 2015, 54, 4495
4. Andreja Abina, Uros Puc, Anton Jeglic, Jana Prah, Rimvydas Venckevicius, Irmantas Kasalynas, Gintaras Valusis, Aleksander Zidansek
"Qualitative and quantitative analysis of calcium-based microfillers using terahertz spectroscopy and imaging" Talanta 2015, 143, 169-177
5. Abina et al, "Applications of Terahertz Spectroscopy in the Field of Construction and Building Materials", Applied Spectroscopy Reviews 50, 279 (2015)
6. G. Montemezzani, M. Alonzo, V. Coda, M. Jazbinsek, P. Gunter, "Running electric field gratings for detection of coherent radiation", J. Opt. Soc. Am. B 32, 1078-1083 (2015)
7. C. Vicario, B. Monoszlai, M. Jazbinsek, O.P. Kwon, C.P. Hauri, "Intense, carrier frequency and bandwidth tunable quasi single-cycle pulses from an organic emitter covering the Terahertz frequency gap", Nature Scientific Reports 5, 14394 (2015)
8. Majkic, M. Zgonik, A. Petelin, M. Jazbinsek, B. Ruiz, C. Medrano, P. Gunter, "Terahertz source at 9.4 THz based on a dual wavelength infrared laser and quasi-phase matching in organic crystals OH1" Appl. Phys. Lett. 105, 141115 (2014)
9. Monoszlai, C. Vicario, M. Jazbinsek, C. P. Hauri, "High-energy terahertz pulses from organic crystals: DAST and DSTMS pumped at Ti:sapphire wavelength" Optics Letters 38, 5106-5109 (2013)

List of publications: THz generation and Detection using THz generators: DAST, DSTMS & OH1

DSTMS

1. Coherent two dimensional-terahertz-terahertz-Raman spectroscopy I-A. Finneran, R. Welsch, M.A. Allodi, T.F. Miller III, and Geoffrey Blake PNAS Vol. 113, 6857-6861, doi:10.1073/PNAS.1605631113; 2016
2. High-performing nonlinear visualization of terahertz radiation on a silicon charge-coupled device,
3. M. Shalaby, C. Vicario, and C.P. Hauri Nat. Commun. 6, 8439 (2015)
4. Ultra-broadband terahertz pulses generated in the organic crystal DSTMS
C. Somma, G. Folpini, J. Gupta, K. Reimann, M. Woerner, and T. Elsaesser Optics Letter 40, 3404, 2015
5. M. Shalaby, C.P. Hauri
Demonstration of a low-frequency three-dimensional terahertz bullet with extreme brightness Nature Comm.2015, 6, No. 5976, doi:10.1038/ncomms 6976
6. C. Vicario, M. Jazbinsek, A.V. Ovchinnikov, O.V. Chefonov, S.I. Ashitkov, M.B. Agranat, and C.P. Hauri
High efficiency THz generation in DSTMS, DAST and OH1 pumped by Cr:forsterite laser Opt. Express, 2015, 23, 4573-4580.
7. M. Shalaby, C.P. Hauri
Terahertz brightness at the extreme: demonstration of 5 GV/m low frequency k terahertz bullet <http://arxiv.org/pdf/1407.1656v1.pdf> (2014)
8. C. Vicario B. Monoszlai, C. P. Hauri
GV/m Single-Cycle Terahertz Fields from a Laser-Driven Large-Size Partitioned Organic Crystal
Phys. Rev. Lett., 2014, 112, 213901
9. C. Vicario, B. Monoszlai, B. ruiz, M. Jazbinsek, C. Medrano and C.P. Hauri Terahertz emission in organic crystals pumped by conventional laser wavelength
SPIE OPTO 89850, 89850C (2014)
10. B. Monoszlai, C. Vicario, M. Jazbinsek, C.P. Hauri.
High-energy terahertz pulses from organic crystals: DAST and DSTMS pumped at Ti:sapphire wavelength
Opt. Lett., 38, No. 23, 5106, 2013
11. Ruchert, C.; Vicario, C; Hauri, C.P.

Spatiotemporal Focusing Dynamics of Intense Supercontinuum THz Pulses
Phys. Rev. Lett., 2013, 110, 123902

12. Vicario, C.; Ruchert, C.; Hauri, C.P.
High field broadband THz generation in organic materials
J. Mod. Opt., 2013, DOI: 10.1080/09500340.2013.800242
13. Stillhart, M.; Schneider, A. & Gunter, P.
Optical properties of 4-N,N-dimethylamino-4'-N'-methyl-stilbazolium 2,4,6-trimethylbenzenesulfonate crystals at terahertz frequencies
J. Opt. Soc. Am. B, 2008, 25, 1914-1919
14. Mutter, L.; Brunner, F. D. J.; Yang, Z.; Jazbinsek, M. & Gunter, P.
Linear and nonlinear optical properties of the organic crystal DSTMS
J. Opt. Soc. Am. B, 2007, 24, 2556-2561
15. Yang, Z.; Mutter, L.; Stillhart, M.; Ruiz, B.; Aravazhi, S.; Jazbinsek, M.; Schneider, A.; Gramlich, V. & Gunter, P.
Large-size bulk and thin-film stilbazolium-salt single crystals for nonlinear optics and THz generation Adv. Funct. Mater., 2007, 17, 2018-2023

DAST

16. C. Vicario, B. Monszlai, G. Arisholm and C.P. Hauri
Generation of 1.5-octave intense infrared pulses by nonlinear interactions in DAST crystal
J. Opt. 17 094005 (2015).
17. C. Vicario, M. Jazbinsek, A.V. Ovchinnikov, O.V. Chefonov, S.I. Ashitkov, M.B. Agranat, and C.P. Hauri
High efficiency THz generation in DSTMS, DAST and OH1 pumped by Cr:forsterite laser
Opt. Express, 2015, 23, 4573-4580.
18. Hauri, C. P.; Ruchert, C.; Vicario, C. & Ardana, F.
Strong-field single-cycle THz pulses generated in an organic crystal
Appl. Phys. Lett., 2011, 99, 161116
19. Hauri, C. P.; Ruchert, C.; Vicario, C. & Ardana, F.
Laser driven generation of intense single-cycle THz field
Proc. SPIE, 2012, 8261, 82610Z
20. Martin, M.; Mangeney, J.; Crozat, P. & Mounaix, P.
Comparison of GaAs and DAST electro-optic crystals for THz time domain spectroscopy using 1.55 μm fiber laser pulses
Terahertz Technology And Applications Iv, 2011, 7938, 793807
21. Cunningham, P. D. & Hayden, L. M.
Optical properties of DAST in the THz range
Opt. Express, 2010, 18, 23620-23625
22. Martin, M.; Mangeney, J.; Crozat, P. & Mounaix, P.
Optical phase detection in a 4-N,N-dimethylamino-4'-N'-methyl-stilbazolium tosylate crystal for terahertz time domain spectroscopy system at 1.55 μm wavelength
Appl. Phys. Lett., 2010, 97, 111112
23. Liu, J. J.; Schmutz, H. & Merkt, F.
Generation of widely tunable Fourier-transform-limited terahertz pulses using narrowband near-infrared laser radiation
J. Mol. Spectrosc., 2009, 256, 111-118
24. Liu, J. J. & Merkt, F.
Generation of tunable Fourier-transform-limited terahertz pulses in 4-N, N-dimethylamino-4'-N'-methyl stilbazolium tosylate crystals
Appl. Phys. Lett., 2008, 93, 131105
25. Jazbinsek, M.; Mutter, L. & Gunter, P.
Photonic Applications With the Organic Nonlinear Optical Crystal DAST
IEEE J. Sel. Top. Quantum Electron., 2008, 14, 1298-1311
26. McLaughlin, C. V.; Hayden, L. M.; Polishak, B.; Huang, S.; Luo, J. D.; Kim, T. D. & Jen, A. K. Y.
Wideband 15 THz response using organic electro-optic polymer emitter-sensor pairs at telecommunication wavelengths
Appl. Phys. Lett., 2008, 92, 151107
27. Schneider, A.; Brunner, F. D. J. & Gunter, P.
Determination of the refractive index over a wide wavelength range through time-delay measurements of femtosecond pulses
Opt. Commun., 2007, 275, 354-358
28. Zheng, X. M.; McLaughlin, C. V.; Cunningham, P. & Hayden, L. M.
Organic broadband terahertz sources and sensors
J. Nanoelectron. Optoelectron., 2007, 2, 58-76

29. Schneider, A.; Stillhart, M.; Yang, Z.; Brunner, F. & Gunter, P.
Improved emission and coherent detection of few-cycle terahertz transients using laser pulses at 1.5 μm - art. no. 658211
Nonlinear Optics and Applications II, 2007, 6582, 658211
30. Schneider, A.; Neis, M.; Stillhart, M.; Ruiz, B.; Khan, R. U. A. & Gunter, P.
Generation of terahertz pulses through optical rectification in organic DAST crystals: theory and experiment
J. Opt. Soc. Am. B, 2006, 23, 1822-1835
31. Schneider, A.; Stillhart, M. & Gunter, P.
High efficiency generation and detection of terahertz pulses using laser pulses at telecommunication wavelengths
Opt. Express, 2006, 14, 5376-5384
32. Schneider, A. & Gunter, P.
Spectrum of terahertz pulses from organic DAST crystals
Ferroelectrics, 2005, 318, 83-88
33. Zheng, X.; Wu, S.; Sobolewski, R.; Adam, R.; Mikulics, M.; Kordos, P. & Siegel, M.
Electro-optic sampling system with a single-crystal 4-N,N-dimethylamino-4'-N'-methyl-4-stilbazolium tosylate sensor
Appl. Phys. Lett., 2003, 82, 2383-2385
34. Laveant, P.; Medrano, C.; Ruiz, B. & Gunter, P.
Rainbow photonics - Growth of nonlinear optical DAST crystals
Chimia, 2003, 57, 349-35

OH1

35. Two-frequency pulsed YLiF₄ cut of the principal axes and THz generation Alain Brenier
Opt. Letters, 2015, 40, 4496-4499.
36. M. Shalaby, C.P. Hauri
Demonstration of a low-frequency three-dimensional terahertz bullet with extreme brightness
Nature Comm. 2015, 6, No. 5976, doi:10.1038/ncomms 6976.
37. C. Vicario, M. Jazbinsek, A.V. Ovchinnikov, O.V. Chefonov, S.I. Ashitkov, M.B. Agranat, and C.P. Hauri
High efficiency THz generation in DSTMS, DAST and OH1 pumped by Cr:forsterite laser
Opt. Express, 2015, 23, 4573-4580.
38. M. Shalaby, C.P. Hauri
Terahertz brightness at the extreme: demonstration of 5 GV/m low frequency k terahertz bullet
<http://arxiv.org/pdf/1407.1656v1.pdf> (2014)
39. A. Majkic, M. Zgonik, A. Petelin, M. Jazbinsek, B. Ruiz, C. Medrano, and P. Gunter.
Terahertz source at 9.4 THz based on a dual-wavelength infrared laser and quasi-phase matching in organic crystals OH1
Appl. Phys. Lett., 2014, 105, 141115
40. A.G. Stepanov, C. Ruchert, J. Levallois, C. Erny and C.P. Hauri
Generation of broadband THz pulses in organic crystal OH1 at room temperature and 10 K
Opt. Mat. Express 4, 870 (2014)
41. Ruchert, C.; Vicario, C.; Hauri, C. P.
Scaling submillimeter single-cycle transients toward megavolts per centimeter field strength via optical rectification in the organic crystal OH1
Opt. Letters, 2012, 37, 899-901
42. Hunziker, C.; Kwon, S. J.; Figi, H.; Juvalta, F.; Kwon, O. P.; Jazbinsek, M. & Gunter, P. Configurationally locked, phenolic polyene organic crystal 2-3-(4-hydroxystyryl)-5,5-dimethylcyclohex-2-enylidenemalononitrile: linear and nonlinear optical properties J. Opt. Soc. Am. B, 2008, 25, 1678-1683
43. Brunner, F. D. J.; Kwon, O. P.; Kwon, S. J.; Jazbinsek, M.; Schneider, A. & Gunter, P.
A hydrogen-bonded organic nonlinear optical crystal for high-efficiency terahertz generation and detection
Opt. Express, 2008, 16, 16496-16508
44. Kwon, O. P.; Kwon, S. J.; Jazbinsek, M.; Brunner, F. D. J.; Seo, J. I.; Hunziker, C.; Schneider, A.; Yun, H.; Lee, Y. S. & Gunter, P.
Organic Phenolic Configurationally Locked Polyene Single Crystals for Electro-optic and Terahertz Wave Applications
Adv. Funct. Mater., 2008, 18, 3242-3250